


## LABORATORY JAW CRUSHER

MODEL NO SL -BT -014, MAKE - LABTEK, ORIGIN - INDIA

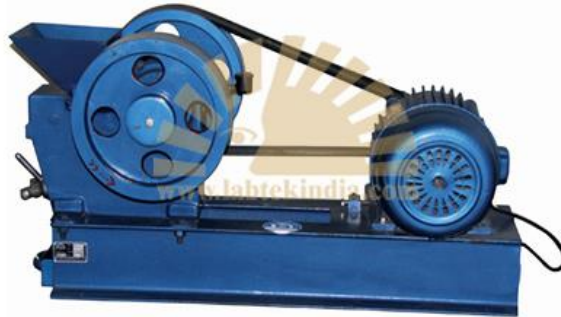

The laboratory jaw Crusher is designed for fast crushing of Aggregates, Ores, Minerals, Coal, Coke, Chemicals and other similar materials. It is compact and of rugged construction for general laboratory of small pilot plant operations.

Two jaws of manganese steel are provided in this laboratory jaw Crusher. The moveable jaw produces two blows for every revolution, thus reducing over sizing to a minimum. A combination of forward and downward strokes with a rocking action exerts pressure on the coarser material, yet permits the finished material to pass through the jaws. A hopper is provided at the top for feeding material. The smooth jaws ensure a uniform product and easy cleaning is possible.